

## Comfort and Convenience

Harmar's SL350AC Stair Lift is engineered with a carriage that is pre-installed on a short section of the durable extruded aluminum track. The backrest has three different settings, and the fold-up footrest adjusts to two different heights for the utmost in custom comfort. The spacious, wide seat features thick padding and attractive upholstery.



### Features & Benefits


- Designed for long-term indoor use
- Constant-pressure rocker controls
- Operates on standard 115 VAC household current
- Seat and footrest fold to save space
- Durable off-white powder coat finish with choices of 3 upholsterys
- 350 lbs. capacity




Optional Grocery Basket



Optional Flip-Up Arms

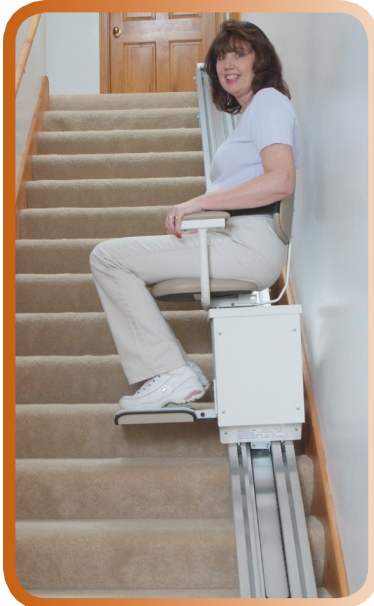

Available Seat Colors:

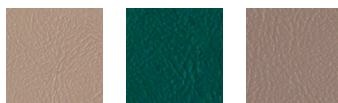

Light Almond Pine Green Mocha  
(Standard)

Specifications	SL350AC
Maximum Weight Capacity:	350 lbs. (160 kg) up to 20 ft.
Overspeed Governor:	No
Track (rail) Type:	Extruded Aluminum
Travel:	16' Standard; 25' Maximum (7.6 m)
Average Number of Return Trips per Charge: (varies with load, length)	N/A
Symmetrical Lift: (Ships and installs non handed)	No
Control in Armrest: (left or right hand operation)	Yes
Need to Open Chassis to Install or Change Sides:	Yes
Lift Mounts to Steps or Wall:	Steps
Minimum Folded Width:	14" (35.5 cm)
Minimum Footrest Height:	4" (10 cm)
Clear Distance Between Armrests:	19" (48.3 cm)
Floor to Seat Height Low / High:	20" - 24" (50.8 cm - 61 cm)
Minimum Wall to Stair Side of Rail:	10" (25.4 cm)
Seat Depth:	14" (35.6 cm)
Backrest Height From Top of Seat:	14" (35.5 cm)
Electrical Requirements:	120 VAC 15A
Operation Power:	AC
Speed:	20 fpm
Incline Range:	30 - 45°
Drive System:	Cable Drive
Warranty:	3-Years on Components
Safety Devices:	Limit Switches Final Limit Switch Slack Cable Brake w/ Switch Footrest Obstruction Switch Swivel Seat Cut-off Switch Constant Pressure Controls Seat Belt