


# PINNACLE STAIR LIFTS

Stepping Up Your Confidence and Independence



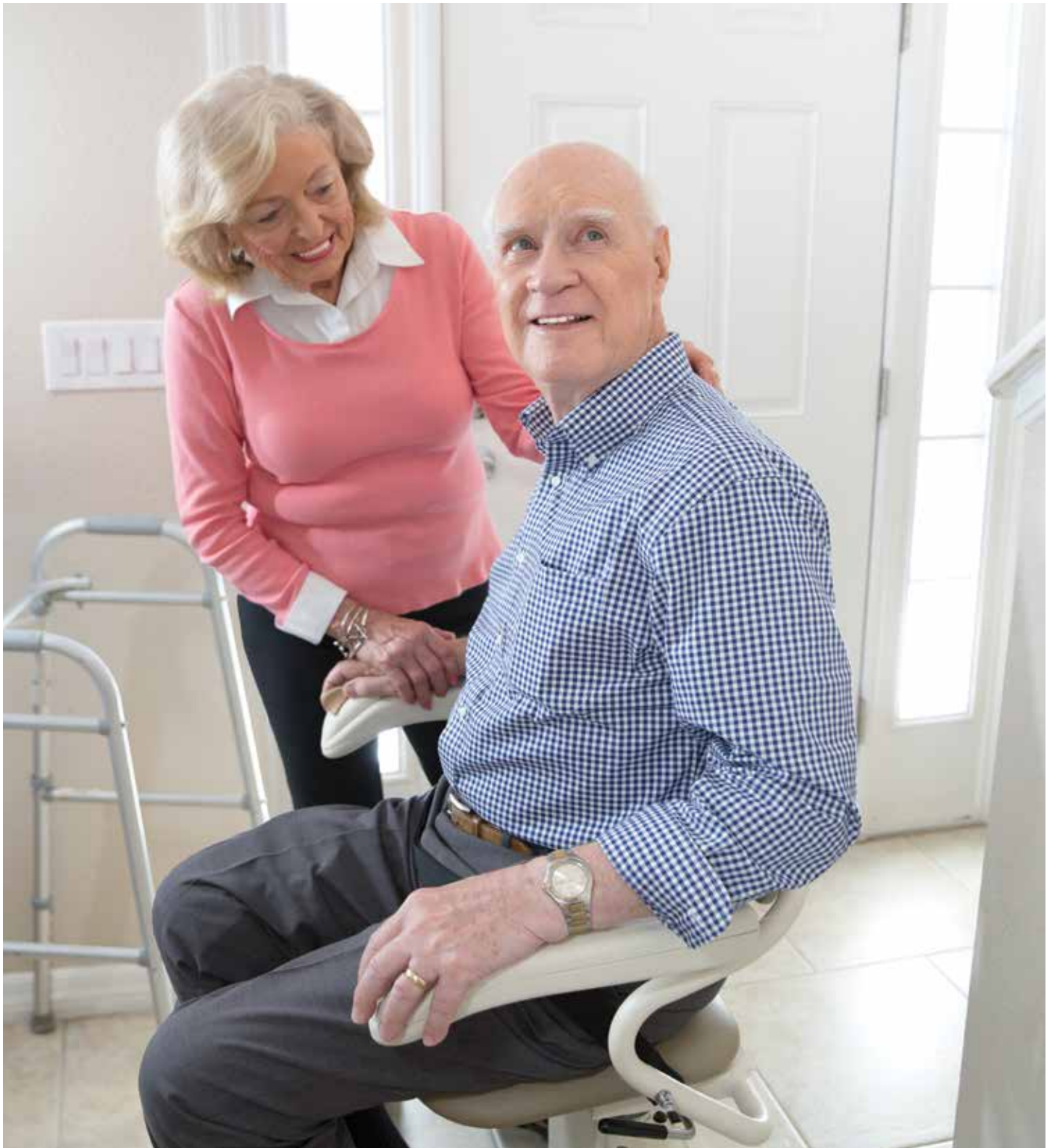

---

Harmar is a leading manufacturer of lifts and is inspired by those who are mobility challenged. Harmar is driven to create mobility and accessibility products that help or enable people to go where they want, when they want, how they want.

Founded in 1998, Harmar has a proud history of growth by evolving our product lines to enrich the lives of our customers. We give people the path to independence to remain in their home, enjoy their families, and be engaged in their communities. Through innovation we are committed to help people of all ages overcome their limited physical mobility at home and away. All Harmar lifts are assembled in the USA, which means quality you can depend on.



# Be Independent. Be Confident. Be Informed.



---

19 million people in the U.S. have disabilities that create problems walking, climbing or descending stairs.

---



---

65% of single family and attached homes in the U.S. have 2 or more levels.

---



---

Every 11 seconds, an older adult is treated in the emergency room for a fall.

---


Whether due to an injury or illness, age or ability, the thought of getting up and down a staircase becomes more challenging. 19 million people in the U.S. have disabilities that create problems walking, climbing or descending stairs... making it more difficult to stay in the home they love. As more individuals plan to age in place, home access is gaining in importance. Safety and security are primary concerns, but so is being able to enjoy your life and family.

Homes are often our foundation of happiness: where good times are had and memories are created. Don't let your staircase become an obstacle to your quality of life.

# What to Look for When Purchasing

There are many stair lifts to choose from in the marketplace. What makes the Pinnacle different?

- Low Maintenance
- Compact
- Helical Worm Gear and Nylon Polymer Gear Rack for a smooth ride and greater stability
- Folding Rail Design
- Efficient motor that gives extended performance during power loss

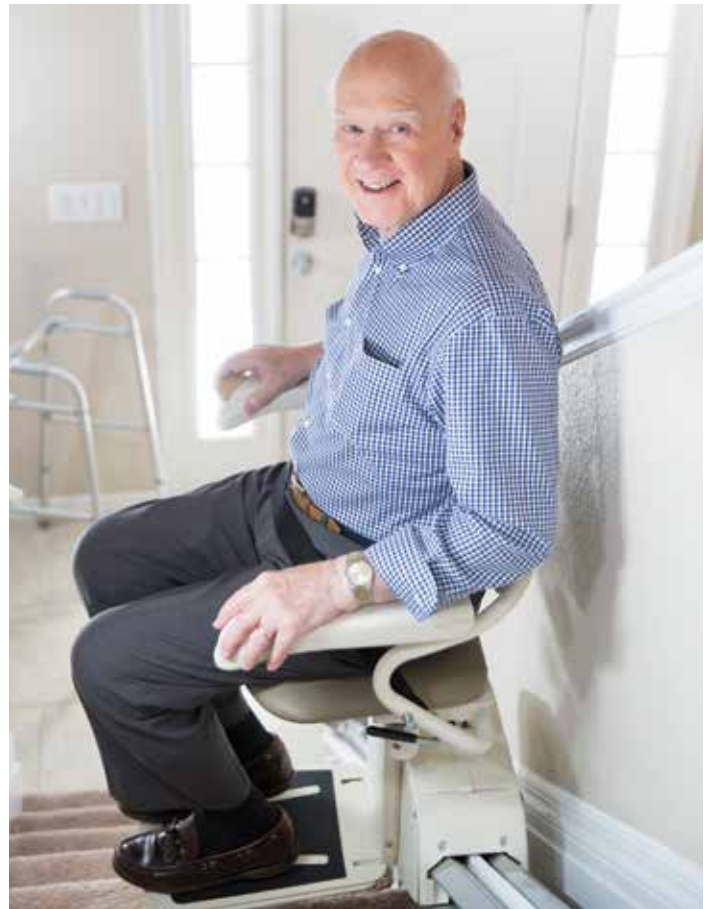

## ZERO MAINTENANCE

For most people that require a stair lift, performing regular maintenance on their lift is a challenge. Harmar's Pinnacle is the only stair lift on the market that does not require lubrication due to its nylon polymer gear rack and worm gear. No mess, no fuss, and no expensive service contracts.



## RAIL SIZE & MOUNTING

Harmar's Pinnacle rail is best in class and sits just 3.5" - 5.5" high. The most compact stair lift on the market, with a folded width of just at 10.5", mounting as close as 0.5" from the wall and allowing for others to easily pass by.



# The Pinnacle Advantage



INDEPENDENT



CONFIDENT



INFORMED

The Pinnacle stair lift gives you the freedom to go up and down stairs with confidence. It gives Independence. It gives Peace of Mind. The Pinnacle isn't just a stair lift, it's a quality of life. For those who own it and their families.



*Every year, nearly 100,000 people make the decision to invest in a stair lift. Our family was one of them. We wanted the same easy access for our son that other families have without thinking of it.*

## Standard Features



Safety obstruction sensors.



Parts warranty



Swivel seat for a safe exit



Wide footrests for stability and extra comfort.\*


*\*Wider foot rest vary slightly in style on each Pinnacle model. See model detail for more information.*

# Peace of Mind Benefits Only with Pinnacle



Many stair lifts on the market have 265 lbs - 300 lbs capacity. While that may be enough for some, there may be need for greater range. More than one family member may need to use the lift.

If these concerns are weighing on you, consider Harmar's SL600 and SL600HD. These lifts boast a 350 lb or 600 lb weight capacity for peace of mind knowing your lift is capable to handle your needs.

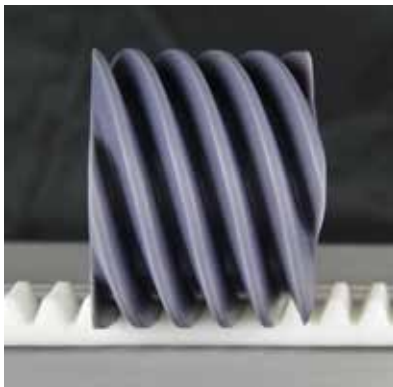

Oversized Helical Worm Gear with six points of contact for an energy efficient low friction design. That means a smooth and stable ride.



Nylon Polymer Gear Rack with steel reinforcement. Lubrication free. That means durability and no track maintenance.



The Pinnacle is the world's most compact stair lift, folding closer to the wall. That means others can pass by easily and more room for you.

# Pinnacle Folding Rail Eliminates Trip Hazards

The Pinnacle Folding Rail deploys automatically to keep walkways clear. This exclusive Harmar design option eliminates the need to stop and push a button or leave the seat. And its simple mechanical design makes it reliable and minimizes system failures.



All Pinnacles come with two remotes so you may operate your chair without sitting on it, or a second user can call the chair if required.

# Harmar Pinnacle Stair Lift Specifications

PINNACLE MODEL	SL600	SL600HD	SL300
<b>Weight Capacity</b>	350 lb	600 lb	300 lb
<b>Drive System</b>	Nylon Polymer Helical Worm Gear	Nylon Polymer Helical Worm Gear	Nylon Polymer Helical Worm Gear
<b>Speed</b>	20 FPM	14 FPM	20 FPM
<b>Seat Backrest</b>	17.5"w x 17"h	17.5"w x 17"h	16.5"w x 9.5"h
<b>Seat Bottom</b>	20"w x 15.5"d	24"w x 16"d	19"w x 15"d
<b>Folded Width</b>	10.5"	17"	11.0"
<b>Folded Width with Folding Rail</b>	13.6"	N/A	13.6"
<b>Control Circuit</b>	24V DC	24V DC	24V DC
<b>Incline Limits</b>	32° - 52°	32° - 52°	32° - 45°
<b>Warranty</b>	3-Year Parts 10-Year Gear Rack 1-Year Battery	3-Year Parts 10-Year Gear Rack 1-Year Battery	2-Year Parts 10-Year Gear Rack 1-Year Battery

## Other Lifts Offered by Harmar



Curved Stair Lifts



SL350 / SL3500D



IL500